









Installation Guide
GCX[®] LIGHT DUTY ROLL STAND
(Refer to qualified personnel)

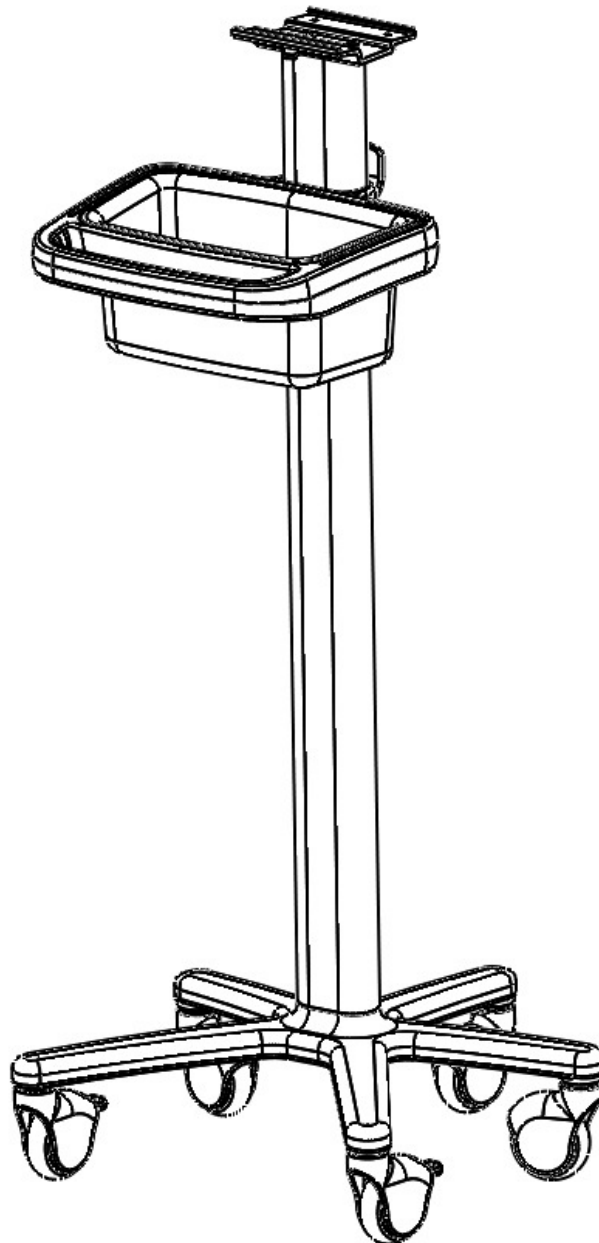

Parts Reference


The following parts and hardware are included in this installation kit:

Item #	Description	Qty.
1	LDRS Column	1
2	LDRS Plunger Plate	1
3	LDRS Base	1
4	LDRS Handle	1
5	LDRS Bin	1
6	M6, 25mm Length, Socket Head Cap Screw	4
7	M6 Split Lock Washer	4
8	M5, 25mm Length, Thread Forming Screw	2
9	M6, 12mm Length, Button Head Cap Screw, Patch	4

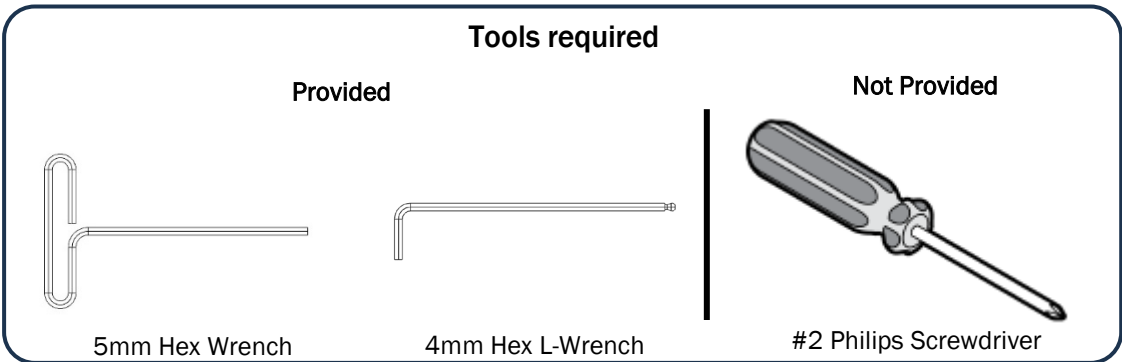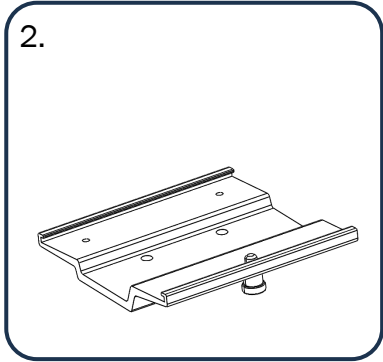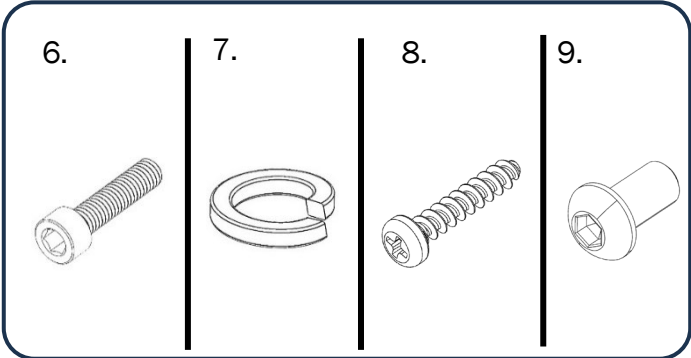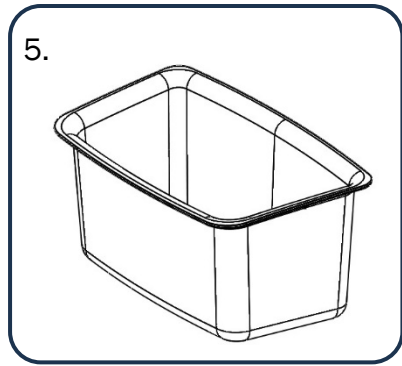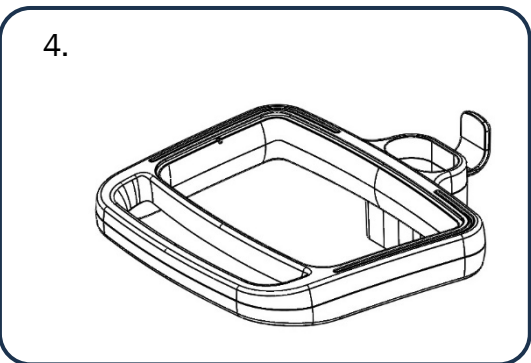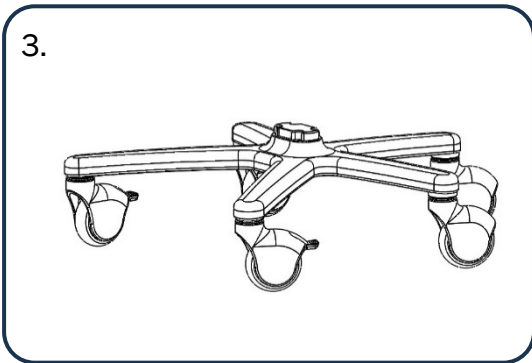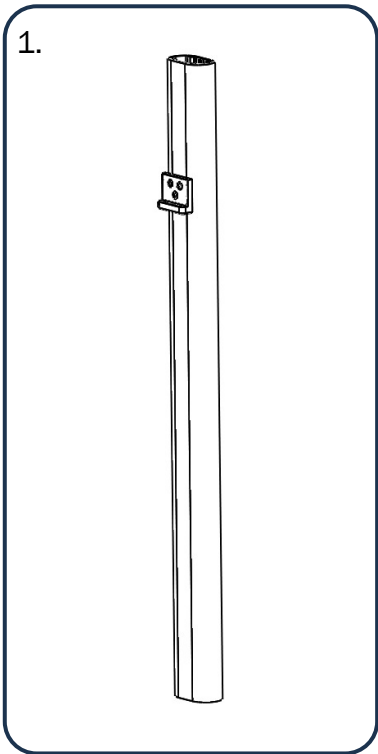


Assembling Column to Base

[Step description]

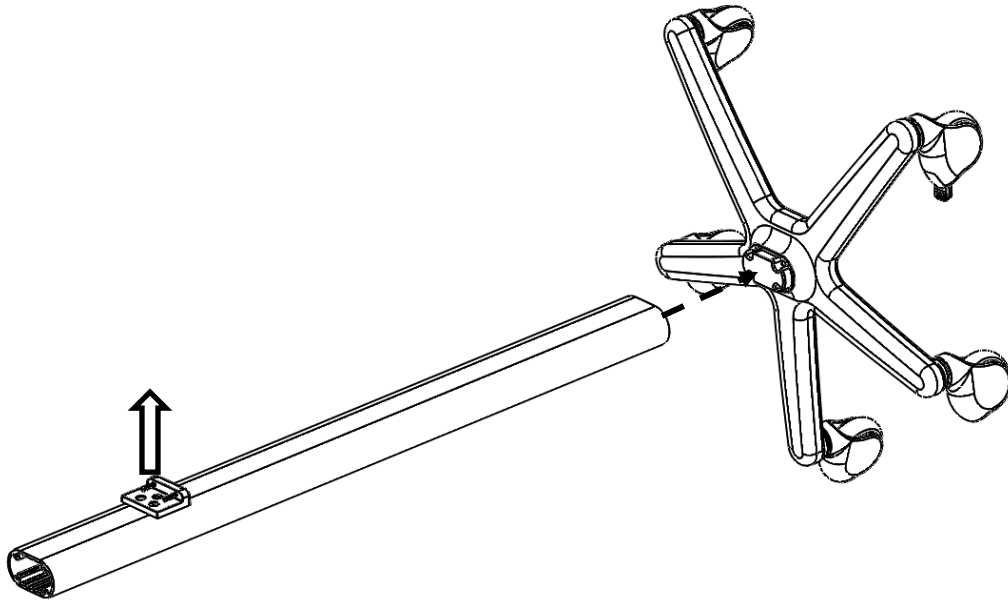


5mm Hex Wrench

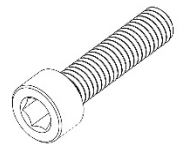


(4) M6 x 1, 25mm Length,
Socket Head Cap Screw




(4) M6 Split Lock Washer

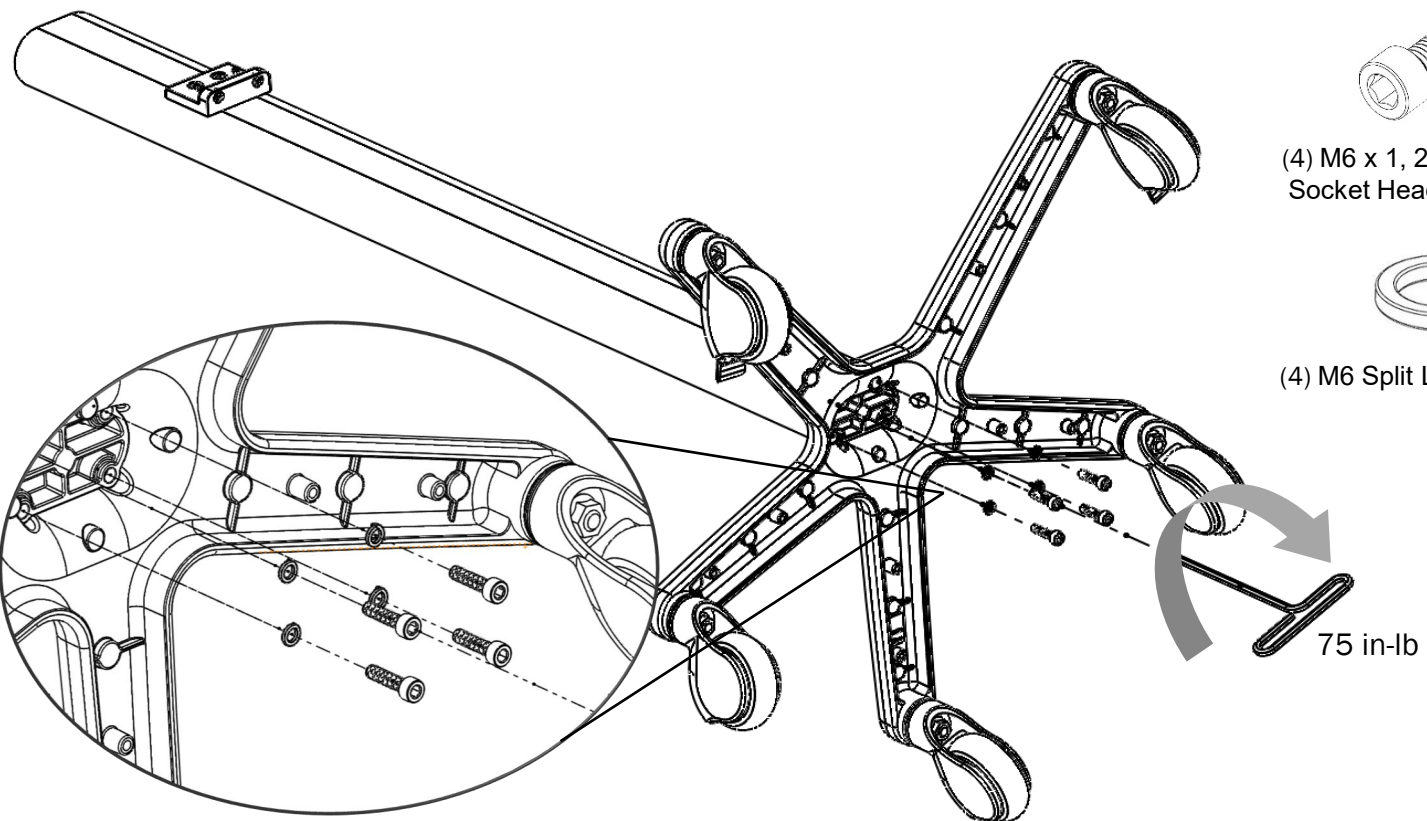

Assembling Handle to Column

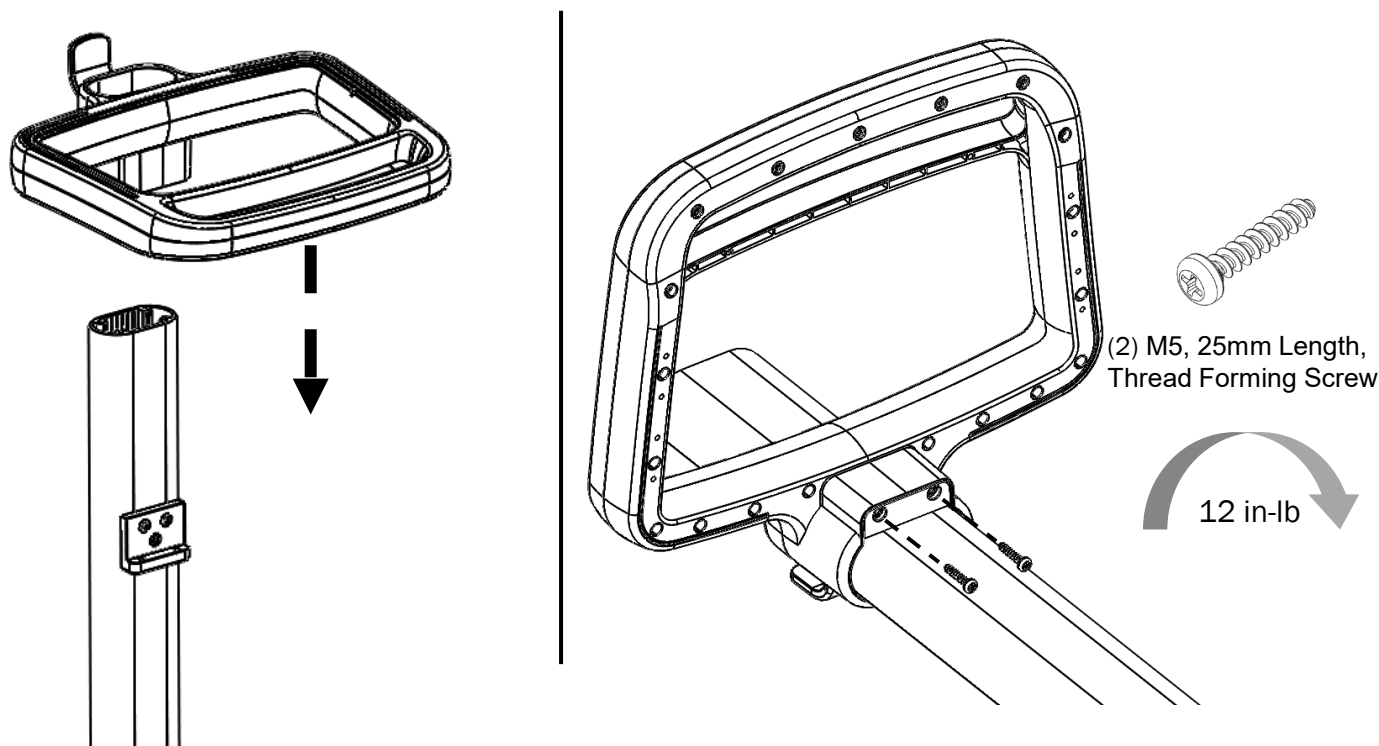


Assembling Plunger Plate to Column

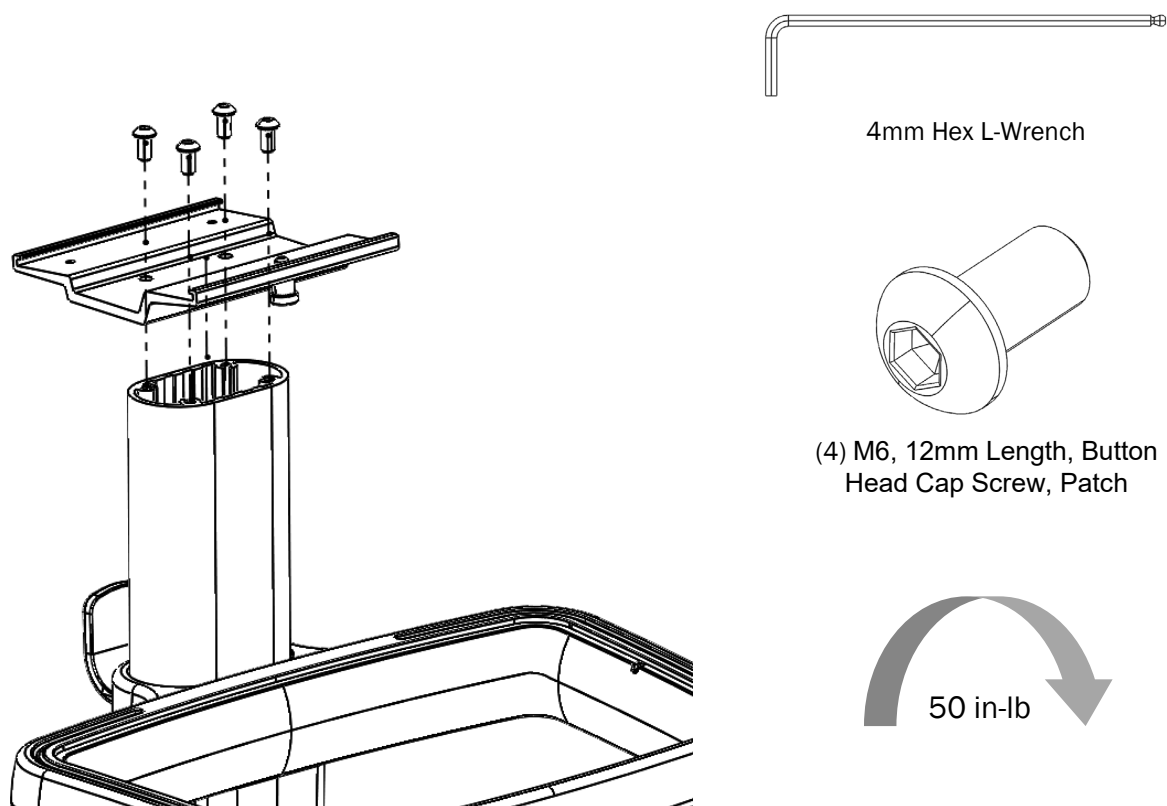


Max Load (W)

6.8 kg 15 lb

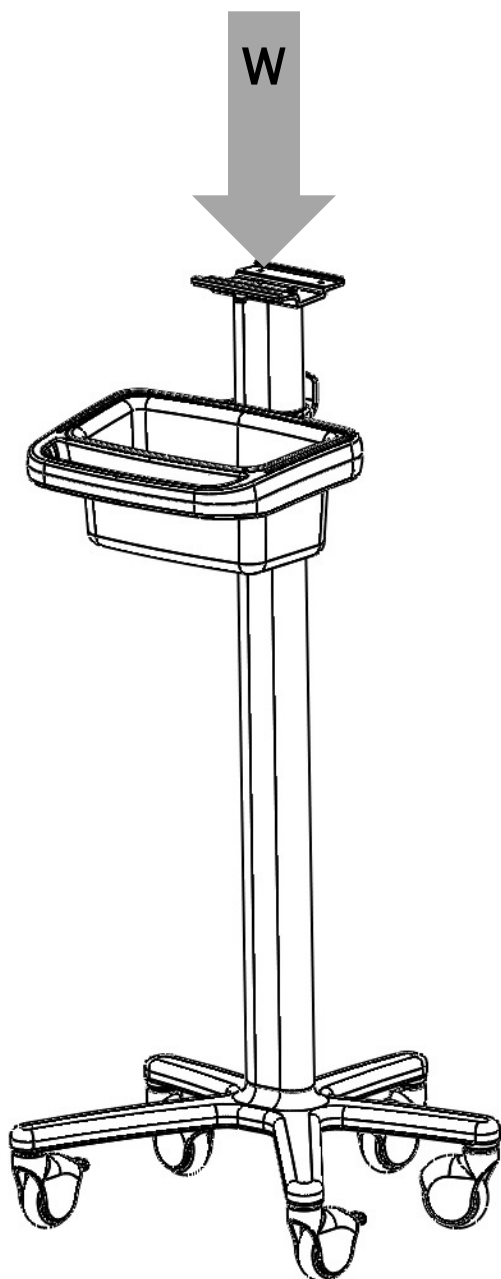


Max Load (W)

2.3 kg 5 lb

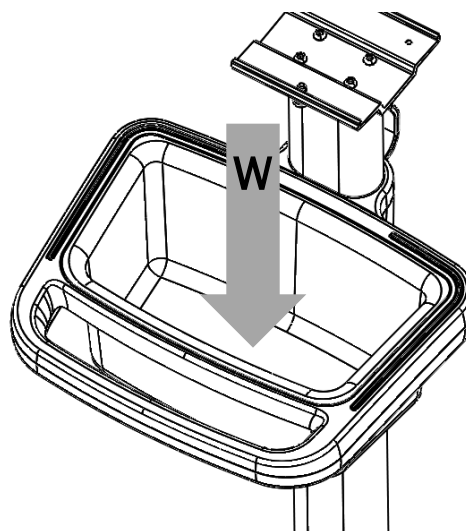

Max Load (W)

2.3 kg 5 lb

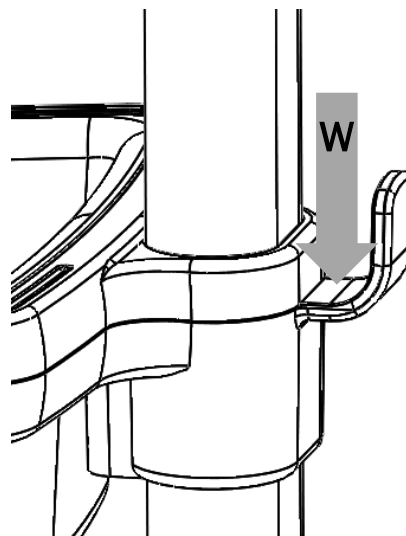

Routine Maintenance

Periodically inspect all fasteners associated with the mounting assembly. Tighten or adjust as necessary for optimal operation and safety.

Cleaning the Mount

The Mount may be cleaned with most mild, non-abrasive solutions commonly used in the hospital environment (e.g. diluted bleach, ammonia, or alcohol solutions).

The surface finish will be permanently damaged by strong chemicals and solvents such as acetone or trichloroethylene.

Steel wool or other abrasive material should never be used.

Damage caused by the use of unapproved substances or processes will not be warranted. We recommend testing of any cleaning solution on a small area of the mounting assembly that is not visible to verify compatibility.

Never submerge or allow liquids to enter the mounting assembly. Wipe any cleaning agents off the mounting assembly immediately using a water-dampened cloth. Dry the assembly thoroughly after cleaning.



CAUTION: GCX makes no claims regarding the efficacy of the listed chemicals or processes as a means for controlling infection. Consult your hospital's infection control officer or epidemiologist. To clean or sterilize mounted devices or accessory equipment, refer to the specific instructions delivered with those products.