

EVO™

Zero maintenance. No service contract required.
Industry-leading 10-year warranty.





The EHR cart has evolved! The result of countless customer inputs and continuous design innovation, EVO is changing the way hospitals think about point-of-care equipment investment – and rewriting industry standards for nurse comfort, trouble-free ownership and proven sustainable value.

- **Built to last from aircraft-grade metal**
EVO stands up to whatever your workflows require.
- **Starts at just 72 lbs. with full-shift power**
Perfectly balanced to roll with ease; relieves nurses' backs, shoulders and arms.
- **Instant, silent, one-handed height adjustment**
Comfort Glide Lift Technology™ lets anyone raise or lower the worksurface with just one hand.
- **Superior infection control; UV-ready**
Smooth, high-gloss antimicrobial finish on all metal surfaces withstands medical-grade chemical and UV disinfection.
- **Zero maintenance, by design**
No maintenance required beyond regular cleaning
- **No service contract needed**
Fewer than 1% of EVOs are sold with service contracts or extended warranties.
- **Industry-best 10-year limited warranty**
EVO is engineered for 10+ years between refreshes.



Comfort Glide Lift Technology™: The Innovation Behind the Evolution

EVO's Comfort Glide creates optimal leverage for smooth, one-handed height adjustment. It trims away the bulk that weighs down straight-post carts, and distributes weight evenly for a cart that's easier for nurses to roll and steer. It eliminates the noise, weight and failure point of powered lift. And it's so reliable, we can back EVO with a 10-year warranty.

EVO is easier on your nurses.

Equipped with full-shift on-board or hot-swap LiFePO₄ power, EVO starts at just 72 pounds. That means EVO could weigh 45% less than carts your nurses are using today.

NEW ClearView Storage: See and feel the difference.

Side-opening drawers let nurses view and access contents without retracting the keyboard or stepping back from the cart. [See page 7.](#)



EVO is easier on your IT staff and budget, too.

Faster deployment, fewer service calls, far less downtime.

Fully assembled, easy to integrate, maintenance free and reliable by design, EVO lets your IT staff spend a lot less time with carts – and a lot more time on mission-critical IT initiatives.

Unmatched investment protection.

You've got a lot more than computers riding on your EHR cart investment. EVO keeps your costs low, your clinical workflows efficient, and your patient satisfaction high.

Choose Your Power

Rechargeable On-Board LiFePO₄ Power

EVO positions the inverter/charger high on the cart away from floor-borne dust and dirt that can clog or damage the system. It's a simple, smart design change that yields big advantages.



Simpler, less-frequent maintenance

Any maintenance can be performed from a comfortable seated or standing position, without hoisting the cart.

LiFePO₄ safety

Unlike other lithium batteries, Jaco's lithium iron phosphate (LiFePO₄) cells are not vulnerable to thermal runaway while charging or fully charged.

Fast charging

Batteries charge to 90% capacity in 2 hours or less.

An extra measure of infection control

Jaco's raised power system electronics do not recirculate floor-borne microbes and pathogens back into the air.

4,000 cycles, guaranteed

Over 95% of Jaco LiFePO₄ batteries 3-6 years old are still in use.

Unmatched reliability

On average, Jaco customers report problems with fewer than .25% of Jaco plug-in power systems per year.

Choose Your Power

Dual-Battery LiFePO₄ Hot-Swap

Jaco's dual-battery hot-swap is a 'true' hot-swap solution: The cart can run for hours on the second battery after the first battery runs down. Compare to single-battery solutions with back-up power that gives nurses minutes or less to swap.

Swapping never interrupts care

Hours-long swap window gives nurses the flexibility to swap on their schedule.

Industry-fastest charge time

Charges to 90% capacity in less than 2 hours.

4,000-cycle battery warranty

Big supply chain savings vs. Li-ion hot-swap - purchase far fewer batteries over time.

Easily monitor battery capacity and life

An LED on the battery charger clearly and quickly indicates when a battery is nearing end-of-life. (Optional software is available.)

Safe for nurses and patients

Jaco's LiFePO₄ battery chemistry is not susceptible to thermal runaway, a cause of fire in other lithium battery types.

Cover your workflows with fewer carts

Customers estimate that Jaco hot-swap has reduced their fleet size as much as 30%.



Locate chargers anywhere

Jaco wall-mount and tabletop (shown) chargers are UL certified; 4" charger depth meets typical fire code requirements for hallway mounting. (Optional on-cart charging available.)

Call (800) 649-2278 for a Hot Swap Device Assessment.

Options & Accessories

New EVO ClearView Storage™



Storage designed specifically for nurse comfort! Compact, side-opening EVO ClearView Storage lets nurses view and access drawer contents with less bending and stretching – and without retracting the keyboard or stepping back from the cart. Nurses can work comfortably and confidently between patient beds, at the med dispenser or in other crowded areas. Choose up to 6 ClearView drawers per cart.

ClearView Cabinet Options

Single Drawer

14.5" w x 9.25" d x 4.5" h



Double Drawer

14.5" w x 9.25" d x 8.25" h



Triple Drawer

14.5" w x 9.25" d x 11.9" h



ClearView Drawer and Bin Options*

Single Drawer

12.5" x 7.5" x 3.5"



1 Large Bin

10" x 7" x 3.2"



2 Long Bins

10.5" x 3" x 3"



3 Mini Bins (NEW)

3.5" x 7.3" x 2.75"



Deep Drawer

12.5" x 7.5" x 7.1"



**All drawer and bin dimensions are interior dimensions.*

Options & Accessories

More EVO Options and Accessories



Wire baskets



Wipes holder



Utility cup




Barcode scanner holders

EVO's new design lets you mount baskets, wipes holders and other accessories closer to the front of the cart, where nurses can reach them easily.

Make it *yours.*

Add custom colors and your logo. Reinforce your brand, identify carts by department, or fit your cart to your care environment.

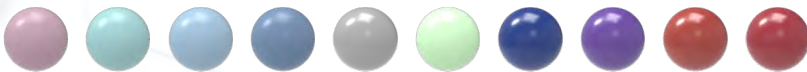

Call (800) 649-2278 for information on all Jaco EVO accessories – or to customize an EVO cart to meet your specific clinical workflow requirements.

Warranties, Service & Support

Expect more. And get more than you expect.



Nobody protects your cart investment like Jaco. We start with EVO's first-ever, industry-leading 10-year warranty. Then we add a 7-year warranty on all Jaco options and accessories, and a 3-year/4,000-cycle warranty on all LiFePO₄ batteries. We even extend our suppliers' part warranties to 3 years. EVO's unmatched warranties and Jaco's track record of proven reliability ensure minimal downtime and maximum return on your investment.

The best service plan is a trouble-free product

Jaco medical computer carts are so reliable, and our warranties so comprehensive, over 99% of our customers do not feel the need to purchase expensive service coverage. Jaco quality can save your hospital hundreds or thousands of dollars on every cart in your fleet.

Jaco Promise Care:

Standard Support with Every Product

- Support from engineers who build our products
- Response in 2 hours or less during business hours (8 a.m. – 5 p.m. ET)
- 24x7 online access to product and service documentation
- Full in-warranty part availability
- Proactive fleet performance reviews
- Product training for your nurses and IT staff



EVO Podium:
The foundation for all EVO models with our industry-best 10-year warranty.

Optional Support and Services

- Factory or onsite cart integration
- Unloading and unpacking service
- Preventative maintenance
- Extended warranties and premium service coverage plans
- Pre-purchased, non-expiring 'blocks' of onsite service time

EVO™



Features & Specifications

Key Features

- Zero-maintenance design
- No service contract needed
- Starting weight just 72 lbs. with full-shift power
- Comfort Glide Lift Technology™ for instant, silent, one-handed height adjustment
- Smooth, high-gloss antimicrobial finish on all metal surfaces withstands medical-grade chemical and UV disinfection
- Choice of on-board LiFePO₄ rechargeable power or dual-battery LiFePO₄ hot-swap power
- Available EVO ClearView Storage™ - side-opening drawer system for improved nurse comfort and convenience
- Durable lightweight metal and steel construction
- Ships assembled and cable-ready
- Manufactured in USA

Warranty

- 10-year limited warranty
- 7-year limited warranty on Jaco manufactured options and accessories
- 3-year limited warranty on electronic components
- 3-year/4,000-cycle limited warranty on all LiFePO₄ batteries

Models

Model EVO-20 LCD

- Worksurface: 20" w x 18" d
- Electronics bay: 16" w x 13" d x 3" h
- Worksurface height range: 33" to 46.5"
- Base: 13" w (front) x 16" d x 8.5" h
- LCD monitor mount:
Left-right pan, 180°;
Tilt up, 90°; Tilt down, 5°;
Maximum LCD size, 24 in.;
Maximum LCD weight, 16 lbs.

(Optional adjustable and heavy-duty LCD mounts available.)

Model EVO-10 Laptop Locking

- Worksurface: 20" w x 16" d
- Screen slot: 16.5" w x 2.75" d
- Electronics bay: 16" w x 13" d x 3" h
- Worksurface height range: 33" to 46.5"
- Base: 13" w (front) x 16" d x 8.5" h

Model EVO PODIUM Podium Cart

Starting weight just 38 lbs.
Contact Jaco for full specifications.

Power Options

Plug-in Rechargeable Power

- Battery chemistry: Lithium iron phosphate (LiFePO₄)
- Output: 120V pure sine wave AC, 180W
- Capacity: 250 Wh or 500 Wh
- Recharge time: 1.75 hours to 80%, 2.5 hours to 100%
- Battery lifespan: 4,000+ recharge cycles per battery
- Certification: UL 60601 (inverter/charger)

Hot-Swap Power

- Battery chemistry: Lithium iron phosphate (LiFePO₄)
- Output: 120V pure sine wave AC, 120W
- Capacity: 430 Wh (with two batteries fully charged)
- Recharge time: ≤ 2.5 hours to 100%
- Battery lifespan: 4,000+ recharge cycles per battery
- Certification: UL 60601 (inverter/charger)



Attend a Live Factory Tele-Demo

Learn the unique way EVO is designed and built, see its features and accessories in action, and get answers to all your questions – in a live private teleconference hosted at Jaco's corporate headquarters.

For more information:

Call **(800) 649-2278**

Visit **www.jacoinc.com/book-teledemo/**

Simple. Solid. Smart.

