



Medical IT Wall Arms and Wall Mounts

A comprehensive portfolio of wall mounted IT workstations designed specifically for patient care

The GCX advantage

50+ years
in hospitals worldwide

A trusted partner at the patient bedside.

Doctors, nurses and clinicians depend on GCX wall arms and mounts for instant, ergonomic access critical health data and information at the point of care. Engineered to outlast and outperform in the most demanding clinical environments, our IT solutions prove themselves every day in ORs, ERs, ICUs and patient rooms across the globe.

Easier to use for nurses and clinicians

Patented ergonomic innovations for effortless adjustment to the most comfortable working position.

More flexible and adaptable

A wide range of options and accessories support multiple clinical applications and workflows.

Safer for everyone

Concealed cables and chemically-resistant components enable superior infection control while eliminating tangled wires and trip hazards.

Less downtime and fewer replacements

Durable construction, extensive testing and industry-best warranties mean fewer interruptions and more value for your investment.

VHM-P Wall Arm

Height-adjustable wall arm for standard and specialty medical IT workstations



True medical-grade arm—load and stability tested to IEC-60601-1 3rd edition standards

16 in. (40.7 cm) 'sit-to-stand' height adjustment plus swivel, tilt and pivot adjustments

Gas-spring counterbalance for effortless height adjustment

Four-bar adjustment mechanism automatically maintains desired viewing angle

Over 80% cable concealment with easy-to-remove/replace cable covers for easy cleaning and maintenance

Aluminum construction with other chemically-resistant materials for long life and superior infection control

Other popular configurations



Slide-out combo keyboard tray/worksurface



Dual monitor (side by side)




Dual monitor (stacked)

VHM-P Arm Accessories



Barcode Scanner Mounts

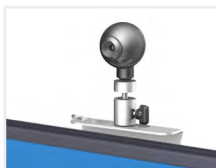

Camera Mounts



Worksurfaces and Trays



Gel Bottle/Probe Cup



USB Charger



Baskets and Bins



EVO Wall Arm

Cost-effective wall arm upgrade

- Comfort Glide Lift Technology™ for effortless height adjustment
- Fits many existing wall channels—upgrade without wiring or wall prep
- Folds to 7.5 in./19.1 cm from wall (8.5 in./21.6 cm with extension arm)
- Optional retractable keyboard (shown) includes spacious worksurface
- Ships fully assembled for fast, easy installation



VHM-25 Variable Height Arm

Space-saving height-adjustable wall arm

- Patented internal mechanism for smooth, counterbalanced adjustment
- Keyboard and screen adjust independently
- Smaller profile—ideal for tight spaces
- Folds close to wall—5 in./12.7 cm with VESA mount, 5.5 in./13.97 cm with folding keyboard



M Series Flush Mount

Versatile flush mount workstation

- Mounts tight to wall but allows tilt and swivel adjustments
- 30 lb./13.6 kg capacity—holds small flat panel display, lightweight keyboard/mouse tray, thin/zero-client PC, accessories
- Multiple mounting options for displays, worksurfaces, peripherals and more
- Mounts to walls, vertical rails, and vertical or horizontal posts



EVO Wall Cabinet

Big-screen in-wall or flush mount workstation

- Holds flat-panel display up to 27 in./68.6 cm—ideal for EHR and medical imaging
- Mounts in wall (projects just 1 in./2.5 cm) or flush on wall (4 in./10.6 cm)
- Height-adjustable keyboard for working comfort
- Two-piece design simplifies wiring and installation

Wall Mount Specifications

	VHM-P	VHM-25	EVO Wall Arm	M Series Flush	EVO Wall Cabinet
Height adjustment range	16 in. (40.6 cm)	12 in. (30.5 cm)	15.75 in. (40 cm)	n/a	n/a
Max. display size	n/a	n/a	24 in. (61 cm)	n/a	27 in. (68.6 cm)
Max. display weight limit	16.6 - 26.6 lbs. (7.5 - 12.1 kg) ¹	3 to 8 lbs. (1.4 - 3.6 kg) or 7 to 20 lbs. (3.2 to 9.1 kg) ^{1,2}	20 lbs. (9.1 kg)	30 lbs. (13.6 kg)	n/a
Dual display options	Side-by-side, stacked	Butterfly, side by side, channel, counter top	n/a	n/a	n/a
Swivel (lateral) at wall	+/-90°	+/-90°	+/-90°	+/-90°	n/a
Display Tilt (+/-)	+/-30°	+/-20°	10°/64°	+/-30°	n/a
Display swivel (pan)	n/a	+/-90°	+/-180°	30°	n/a
Worksurface dimensions	18.3 x 9 in. (46.5 x 22.9 cm) or 20 x 9 in. (50.8 x 22.9 cm)	19.4 x 16 in. (49.3 x 40.6 cm)	20 x 9 in. (50.8 x 22.86 cm)	Several options	19 x 9 in. (48.3 x 22.9 cm)
Arm extensions	8 in. (20.3 cm) 14 in. (35.6 cm)	7 in. (17.8 cm) 14 in. (35.6 cm)	15 in. (38.1 cm)	n/a	n/a
Maximum reach std.	39 in. (99.1 cm)	21.8 in. (55.4 cm)	34 in. (86.4 cm) or 42 in. (106.7 cm)	n/a	n/a
w/ extensions	53 in. (134.6 cm)	29 in. (73.7 cm)	49 in. (124.5 cm) or 57 in. (144.8 cm)	n/a	n/a
Folded distance from wall	9.5 in. (24.1 cm)	5.5 in. (13.97 cm)	7.5 in. (19.1 cm) 8.5 in. (21.6 cm) with extension	n/a	1 in. (2.5 cm) embedded 4 in. (10.6 cm) flush mounted

¹ Limits differ for dual monitor configurations. See gcx.com for details.

² Mounting plate dependent. See gcx.com for details.



sales@gcx.com

3875 Cypress Drive
Petaluma, California 94954
800.228.2555

140 Constitution Blvd.
Franklin, MA 02038
800.649.2278

Dependent on manufacturing location, GCX quality management processes are ISO 13485:2016 certified or ISO 9001:2015 certified. GCX, PolyMount, and PolyQuip are registered trademarks, and VHRC, VHC, VHM, VHRS, PRC, M Series, and Stor-Lock are trademarks, of GCX Corporation. All other brand names, product names, or trademarks belong to their respective holders. © Copyright 2025